


**POWERF2™**





Ideal conditions don't exist—ice or  
incline, gravel or grass, slope or slush.  
When the terrain fights back, you need  
a system that pushes forward.

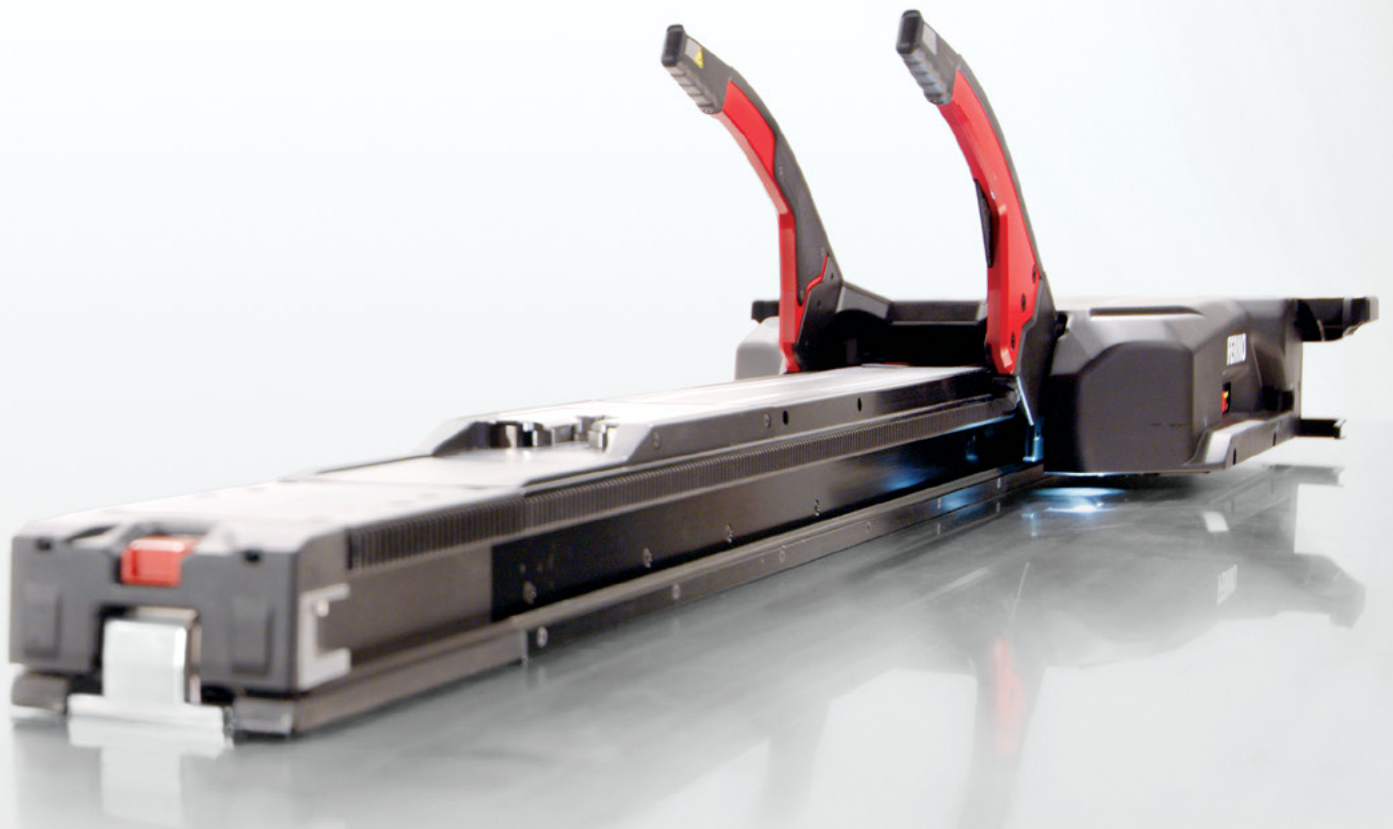




Introducing the all-new FERNO® POWER F2™  
Fastening System with CrossDrive™.

Unlike traditional lift-lower systems, the POWER F2™  
does the work for you, giving emergency personnel  
total control to effortlessly lift, load, and unload  
patients at the touch of a button.

Designed for cross-compatibility, the POWER F2™  
works seamlessly with multiple cot systems, making  
it the most versatile solution available. No need for  
costly fleet upgrades—just more power, control,  
and efficiency when it matters most.

Reduce injuries. Stay in control.  
Be ready for whatever comes next.

LEARN MORE AT  
[POWER-F2.FERNO.COM](https://POWER-F2.FERNO.COM)

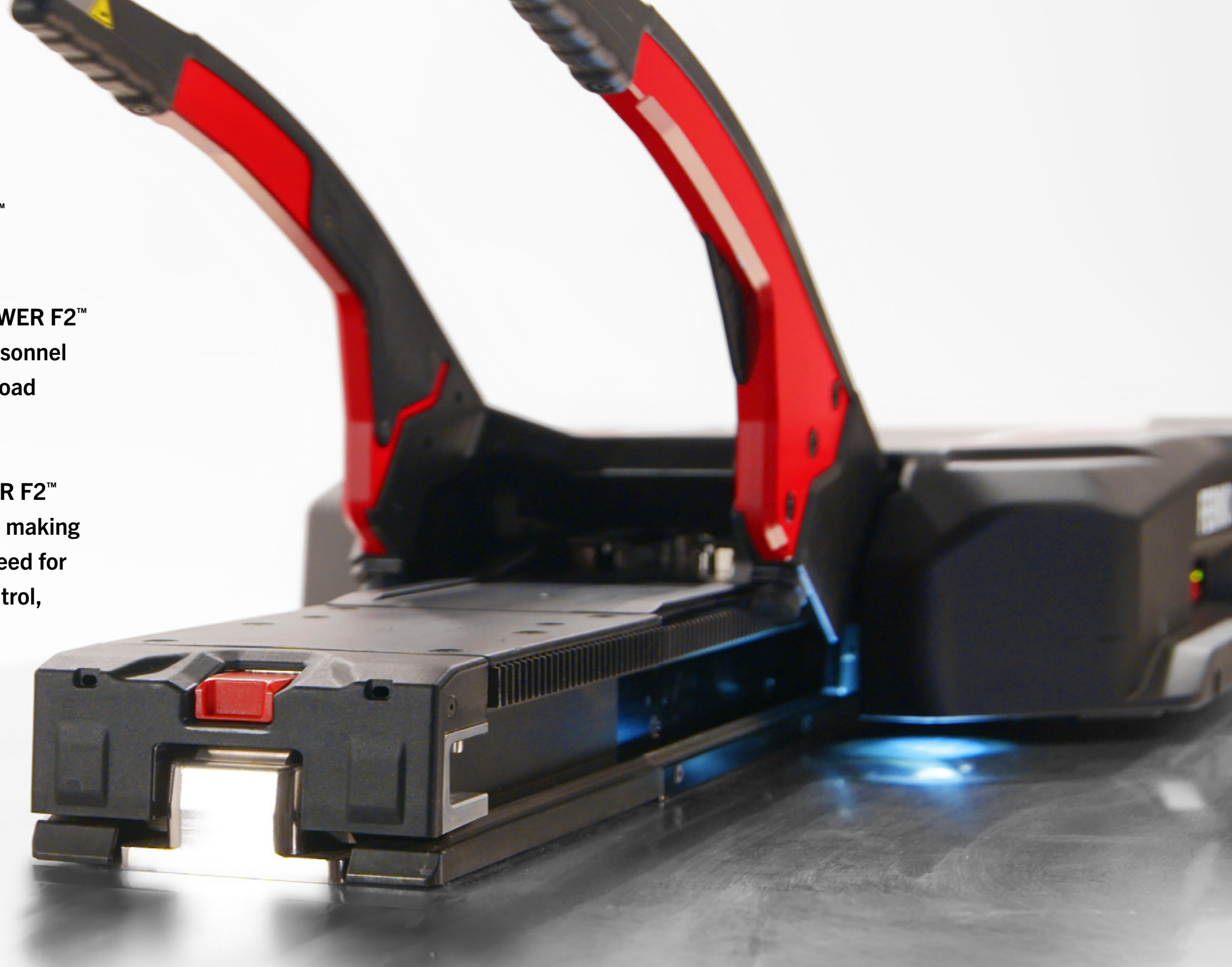





# Powered for efficiency, consistency and longevity.

FERNO's CrossDrive™ technology is transforming the way emergency care is delivered for both caregivers and patients.



- Universally cross-compatible with the majority of manual and powered cots, ensuring seamless integration with existing products
- Power drive mode lifts and loads patients through simple touch controls on the POWER X2™ Cot, FERNO Connect, or POWER F2™ Carriage, minimizing strain and injury risk for paramedics.
- 700 lbs lift capacity eliminates the need for paramedics to bear the weight, ensuring safer and more efficient patient transportation
- Manual mode enables complete loading and unloading of patients for everyday use or in the event of power loss





Eliminating barriers for a whole new powered experience.



- SceneBright LED lights automatically illuminate the ground below the carriage when deployed from the vehicle, ensuring optimal visibility
- Built-in, spare EnduraCharge 36.3 volt LiFePO4 Battery with integrated charging provides a consistent, reliable back-up power for uninterrupted operation
- IP56 rating ensures reliability even in the most demanding environments






# Take control.

The FERNO Connect touchscreen console is the central hub, providing full control and insight into the POWER X2 and POWER F2 systems.

- Action screens in combination with directional LED lights visually guide paramedics during loading and unloading for efficient and safe operation
- Battery charging status for both POWER F2™ and POWER X2™
- Alerts prompt operators and technicians to error states
- Cycle counts for both POWER F2™ and POWER X2™ to help with maintenance and extend equipment life
- Fleet and asset management streamline operations and maximize efficiency
- Flexible and easy mounting options based on vehicle configuration

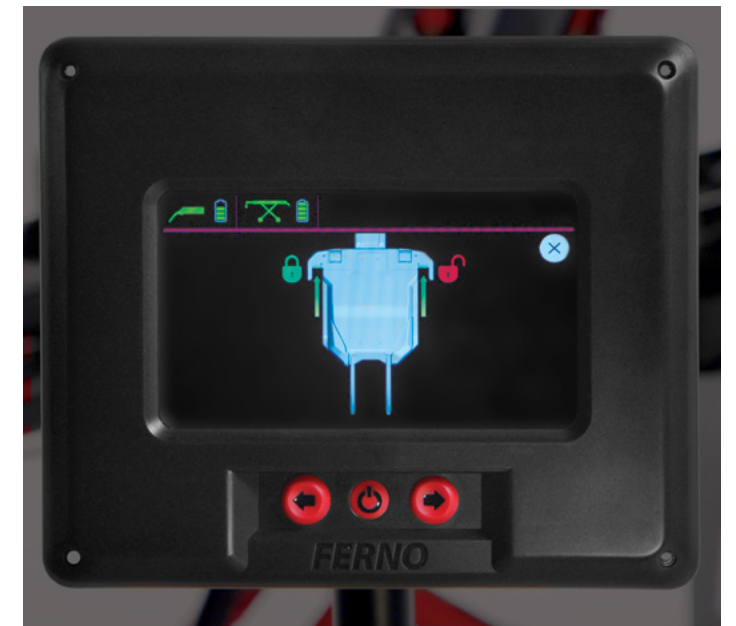

## **POWERX2™**

The new FERNO® POWER X2™ Ambulance Cot is cross-compatible and offers more power, visibility, and maneuverability than today's leading competitive powered cots. It brings stability to unpredictable environments for operators, allowing them to focus on safe patient care and transportation.





## SPECIFICATIONS

OVERALL LENGTH

**88" / 2235 MM**

WIDTH

**24.5" / 622 MM**

WEIGHT

**199.5 LBS / 90.45 KG (TOTAL  
WEIGHT INCLUDING 2 LBS FOR  
FERNO CONNECT)**

**16.5 LBS / 7.5 KG (FLOOR PLATE  
ASSEMBLY)**

**115 LBS / 52.1 KG (FIXED RAIL  
ASSEMBLY)**

**66 LBS / 29.9 KG (CARRIAGE  
ASSEMBLY)**

MAXIMUM LIFT CAPACITY\*

**700 LBS / 318 KG**

RECOMMENDED LOADING HEIGHT

**UP TO 36"**

BATTERY

**36.3V (2.6AH) LITHIUM IRON  
PHOSPHATE (LiFePO4)**

\*MAXIMUM LIFT CAPACITY REPRESENTS PATIENT  
WEIGHT AND ACCESSORIES. FERNO® SAFE WORKING  
LOAD OF 870 LBS (395 KG) REPRESENTS THE SUM  
OF THE COT TOTAL WEIGHT WITH ACCESSORIES  
AND PATIENT

FERNO® RESERVES THE RIGHT TO CHANGE  
SPECIFICATIONS WITHOUT NOTICE.

## COMPLIANCE

**NZS-4535**


**BS EN-1789**

**SAE J3027**

THE FERNO® POWER F2™ IS COMPLIANT WITH  
THE ABOVE STANDARDS WHEN USED WITH:

FERNO® POWER X2™ AMBULANCE COT  
STRYKER® POWER-PRO™ XT AMBULANCE COT  
STRYKER® POWER-PRO™ 2 AMBULANCE COT



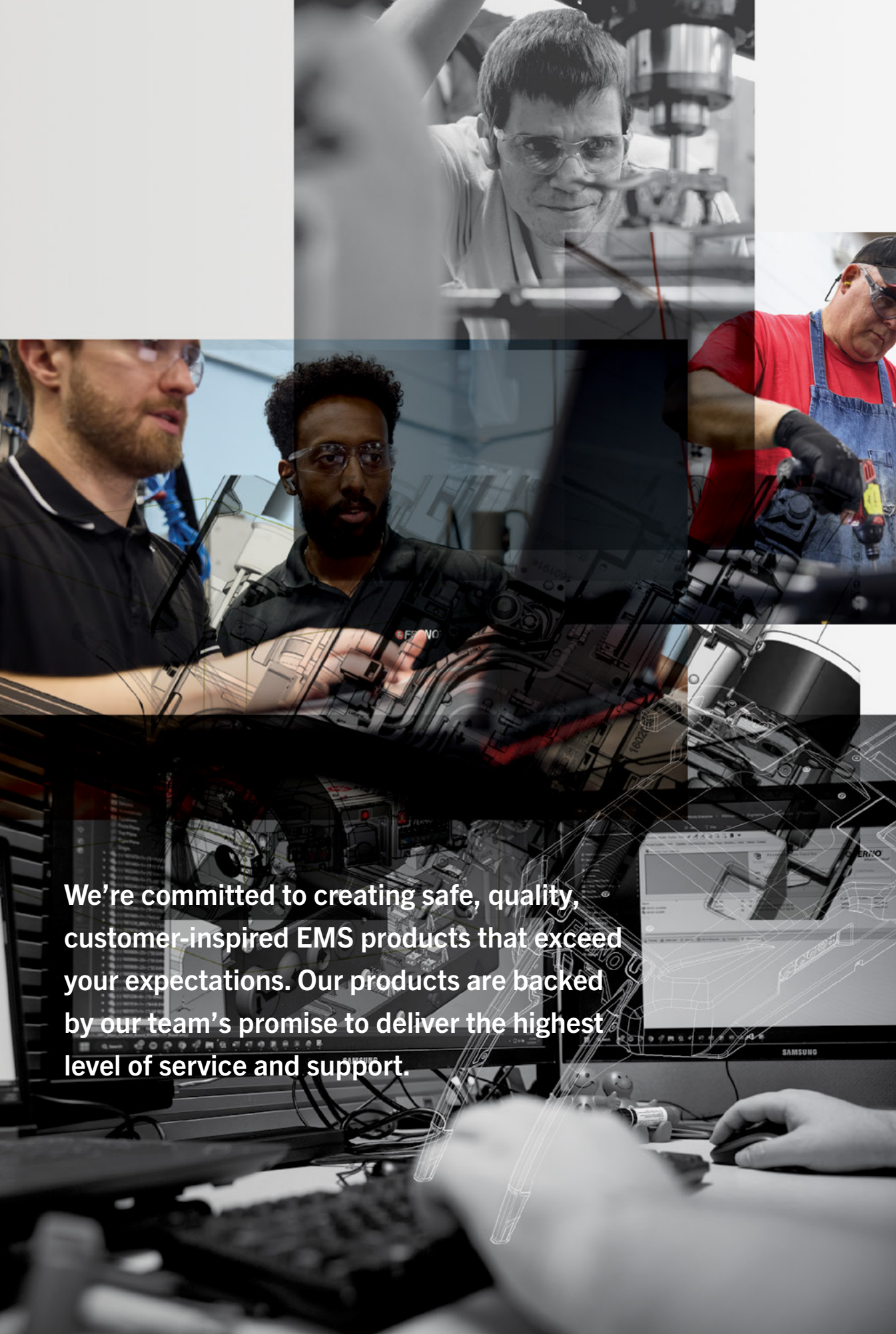

We're committed to creating safe, quality, customer-inspired EMS products that exceed your expectations. Our products are backed by our team's promise to deliver the highest level of service and support.



*"We want to make sure that we're being empathetic and putting ourselves in the user's position... how would they handle a situation and make sure that we do what's best for them."*

NICK  
SENIOR PROJECT ENGINEER, FERNO





FERNO Inc.  
70 Weil Way  
Wilmington, OH 45177

[FERNO.COM](http://ferno.com)

UNITED STATES  
1-877-733-0911  
1-937-382-1191 - Fax  
1-888-388-1349 - Fax Orders  
[fernoorders@ferno.com](mailto:fernoorders@ferno.com)

INTERNATIONAL  
1-937-382-1451  
1-937-382-6569 - Fax  
[intlsales@ferno.com](mailto:intlsales@ferno.com)

SERVICE  
877-733-0911  
[techsup@ferno.com](mailto:techsup@ferno.com)

COT CARE & MAINTENANCE  
[service@ferno.com](mailto:service@ferno.com)

© 2025 FERNO, Inc.  
MAY-25-v3.0