



## Ferno F2 Lift-Off H56 incl mattress & belts




### PRODUCT INFO

Product Nr:	F600331005SN
Length:	197 cm
Width:	58cm
Height:	32-85 cm
Loading:	56 -59 cm
Weight:	46 kg (27,5+18,5)
Max Load:	250 kg
Compliance:	Class 1 MDR / CE
Compliance:	EN1789 & EN1865
Color:	Grey & White

The **Ferno F2 Lift-Off** provides a safe, ergonomic, and reliable patient handling solution. Made of lightweight, durable aluminum, the F2 Lift-Off features an independent leg design that supports the patient during loading and unloading, reducing strain on medics. Automatic locking front and rear legs automatically unfold. The color-coded left- and right-hand control levers are user-friendly. Four three-position, telescoping handles provide ergonomic carrying and lifting options for medics. Adjustable six-position backrest provides support for patient comfort, and optional swing-down side rails provide patients with a feeling of security.

**Ferno F3 Lift-Off** means that the stretcher is two-part. The upper section can be easily detached for carrying the stretcher in stairs or similar situations. The portable stretcher section weighs only 18.5 kilograms

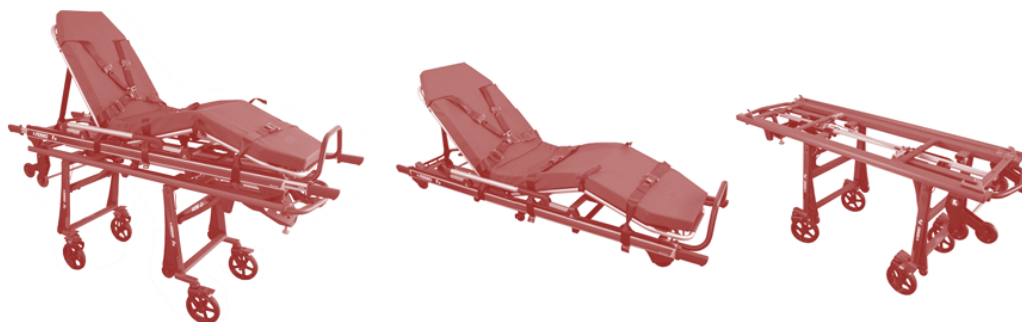

### STANDARD FEATURES

- Excellent maneuverability over any terrain thanks to large 200 mm wheels
- Adjustable multi-position backrest supports a variety of comfort positions
- Fold-downside rails give patients a feeling of security
- Shock frame with contour feature
- Low maintenance
- Includes mattress, patient restraints, sidearms and a support foot
- Telescoping IV-Pole as an option
- Color-coded left-hand and right-hand control levers
- Loading height in Vehicle 56-59 cm

**Option:** Ferno F2 Fastener 10G - F600150003